CHORUS

RUBELLA IgG

REF 81030 (36 tests)

Prodotto da/Manufactured by/Fabricado por:
DIESE Diagnostica Senese
 Via delle Rose 10
 53035 Monteriggioni (Siena) - Italy

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(English)

1. INTENDED USE

QUANTITATIVE DETERMINATION OF IgG-CLASS ANTIBODIES TO RUBELLA VIRUS IN HUMAN SERUM, USING A DISPOSABLE DEVICE APPLIED ON THE CHORUS INSTRUMENTS

2. INTRODUCTION

Rubella virus causes a very common infectious disease in humans, accompanied by a typical diffused rash. The course of the illness is generally short and without complications. However, it is harmful for the fetus if contracted during pregnancy, because the virus is teratogenic and may cause very serious malformations especially during the first month of pregnancy. It is therefore essential to determine the patient's state of immunity before the start of pregnancy, if possible, in order to perform vaccination or, if this is impossible, to monitor the course of the pregnancy to observe whether seroconversion occurs. Several methods are employed for the assay of specific immunoglobulins. Of these methods, the ELISA technique offers the advantage of being completely automated.

3. PRINCIPLE OF THE METHOD

The test is based on the ELISA principle (Enzyme Linked Immunosorbent Assay). The antigen, composed of purified and inactivated Rubella virus, is bound to the solid phase. The specific immunoglobulins are bound to the antigen through incubation with dilute human serum.

After washings to eliminate the proteins which have not reacted, incubation is performed with the conjugate, composed of human IgG monoclonal antibodies conjugated to horse radish peroxidase.

The unbound conjugate is eliminated, and the peroxidase substrate added.

The colour which develops is proportional to the concentration of specific antibodies present in the serum sample.

The disposable devices contain all the reagents to perform the test when applied on the Chorus instruments.


The results are expressed in IU/ml calculated in reference to the 1st International Standard Preparation.

4. KIT COMPOSITION AND REAGENT PREPARATION

The kit is sufficient for 36 tests.

DD DEVICES 6 packages each containing 6 devices

Use: **equilibrate a package at room temperature**, open the package and remove the devices required; replace the others in the bag with the silica gel, expel the air and **seal** by pressing the closure. Store at 2-8°C.

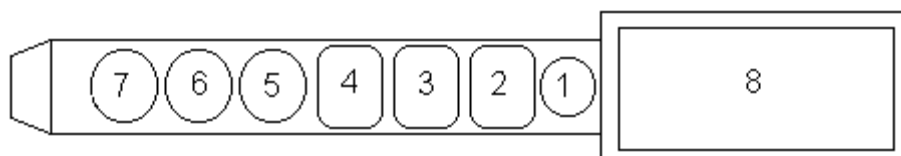Description:

Position 8: Space for application of bar code label

Position 7: empty

Position 6: MICROPLATE WELL coated with Rubella antigen

Position 5: Uncoated MICROPLATE WELL

Position 4: TMB SUBSTRATE

Contents: Tetramethylbenzidine 0.26 mg/mL and H₂O₂ 0.01% stabilized in 0.05 mol/L citrate buffer (pH 3.8)

Position 3: SAMPLE DILUENT

Contents: Proteic solution containing phenol 0.05%, Bronidox 0.02% and an indicator to reveal the presence of the serum.

Position 2: CONJUGATE

Contents: anti-human IgG monoclonal antibodies labelled with horseradish peroxidase, in phosphate buffer containing phenol 0.05% and Bronidox 0.02%.

Position 1: EMPTY WELL in which the operator must place the undiluted serum

CALIBRATOR CALIBRATOR 1 x 0.425 mL

Contents: Diluted human serum, at known antibody concentration, containing Proclin and Gentamycin, **liquid ready for use.**

CONTROL + POSITIVE CONTROL 1 x 0.425 mL

Contents: Diluted human serum, at known antibody concentration, containing Proclin and Gentamycin, **liquid ready for use.**

MATERIALS REQUIRED BUT NOT PROVIDED

- WASHING BUFFER [REF] 83606
- CLEANING SOLUTION 2000 [REF] 83609
- SANITIZING SOLUTION [REF] 83604 – 83608
- CHORUS NEGATIVE CONTROL/SAMPLE DILUENT [REF] 83607
- Chorus / Chorus Trio Instrument
- Distilled or deionised water
- Normal laboratory glassware: cylinders, test-tubes etc.
- Micropipettes for the accurate collection of 50-200 µl solution
- Disposable gloves
- Sodium Hypochlorite solution (5%)
- Containers for collection of potentially infectious materials

5. STORAGE AND STABILITY OF REAGENTS

Reagents must be stored at 2/8°C. In the case of storage at an incorrect temperature the calibration must be repeated and the run validated using the control serum (see point 9, Test validation).

The expiry date is printed on each component and on the kit label.

Reagents have a limited stability after opening and/or preparation:

DEVICES	8 weeks at 2/8°C
POSITIVE CONTROL	8 weeks at 2/8°C
CALIBRATOR	8 weeks at 2/8°C

6. PRECAUTIONS

ONLY FOR DIAGNOSTIC USE IN VITRO.

Attention:

This kit contains materials of human origin which have been tested and gave a negative response by FDA-approved methods for the presence of HbsAg and for anti-HIV-1, anti-HIV-2 and anti-HCV antibodies. As no diagnostic test can offer a complete guarantee regarding the absence of infective agents, all material of human origin must be handled as potentially infectious. All precautions normally adopted in laboratory practice should be followed when handling material of human origin.

Waste disposal: serum samples, calibrators and strips once used must be treated as infectious residuals and eliminated according to law.

Health and Safety Information

1. Do not pipette by mouth. Wear disposable gloves and eye protection while handling specimens. Wash hands thoroughly after placing the devices in the CHORUS instrument.
2. The following reagents contain low concentrations of harmful or irritant substances:
 - a) The conjugate contains phenol
 - b) The substrate is acid
 If any of the reagents come into contact with the skin or eyes, wash the area extensively with water.
3. Neutralized acids and other liquid waste should be decontaminated by adding a sufficient volume of sodium hypochlorite to obtain a final concentration of at least 1%. A 30 minute exposure to 1% sodium hypochlorite may be necessary to ensure effective decontamination.
4. Spillage of potentially infectious materials should be removed immediately with adsorbent paper tissue and the contaminated area swabbed with, for example, 1% sodium hypochlorite before work is continued. Sodium hypochlorite should not be used on acid-containing spills unless the spill area is first wiped dry. Materials used to clean spills, including gloves, should be disposed of as potentially biohazardous waste. Do not autoclave materials containing sodium hypochlorite.

Analytical precautions

Bring the devices to room temperature (18-30°C) before use; use within 60 min.

1. **Discard devices which show the substrate (well 4) blue colored.**
2. Adding the sample into the well verify that it is perfectly distributed on the bottom.
3. Check for the presence of the reagents in the device and that the device is not damaged; do not use devices which are lacking a reagent.
4. The devices are for use with the Chorus instrument; the instructions for use must be carefully followed and the instrument operating manual must be consulted.
5. Check that the Chorus instrument is set up correctly (see Chorus Operating Manual).
6. Do not alter the bar code placed on the handle of the device in order to allow correct reading by the instrument
7. Defective barcodes can be inserted manually in the instrument.
8. Do not expose the devices to strong light or to hypochlorite vapours during storage and use.
9. The use of strongly hemolyzed samples, or samples presenting microbial contamination may all constitute sources of error.
10. Before inserting the devices in the instrument, check that the reaction well does not contain foreign bodies.
11. Pipette the test serum (50 µl) in well 1 of the device (see figure).
12. Do not use the device after the expiry date

7. TYPE AND STORAGE OF SAMPLE

The sample is composed of serum collected in the normal manner from the vein and handled with all precautions dictated by good laboratory practice. The fresh serum may be stored for 4 days at 2/8°C, or frozen for longer periods at -20°C, and can be thawed a maximum of 3 times. Do not keep the samples in auto-defrosting freezers. Defrosted samples must be shaken carefully before use. The quality of the sample can be seriously affected by microbial contamination which leads to erroneous results.

Strongly lipemic, icteric, hemolyzed or contaminated samples should be avoided.

The test cannot be applied to plasma.

8. PROCEDURE

1. Open the package (on the side containing the pressure-closure), remove the number of devices required and seal the rest in the bag after expelling the air.
2. Check the state of the device according to the indications reported in the Analytical Warnings points 1 and 8.
3. Dispense 50 µl of undiluted test serum in well no. 1 of each device; at each change of batch, use a device for the calibrator.
4. Place the devices in the instrument. Perform the calibration (if necessary) and the test as reported in the Chorus Operating Manual.

9. TEST VALIDATION

Use the control serum to check the validity of the results obtained. It should be used as reported in the operating manual. If the instrument signals that the control serum has a value outside the acceptable range, the calibration must be repeated. The previous results will be automatically corrected.

If the result of the control serum continues to be outside the acceptable range, contact the Scientific Support.

Tel: 0039 0577 319554
 Fax: 0039 0577 366605
 email: scientificsupport@diesse.it

10. INTERPRETATION OF THE RESULTS

The Chorus instrument gives a result in IU/ml, calculated on the basis of a lot-dependent graph memorized in the instrument.

The degree of immunity of the test serum can be interpreted as follows:

IMMUNE: when the IgG concentration in the sample is >12 IU/ml

NOT IMMUNE: if the concentration is < 8 IU/ml

DOUBTFUL: when the concentration is in the range 8-12 IU/ml.

In the case of a doubtful result, repeat the test. If the test remains doubtful, collect a new sample after 1-2 weeks.

When the sign > appears on the instrument the IU/mL are over 100 IU/mL; in this case dilute the sample in Sample Diluent/Negative Control **REF** 83607 (e.g. 1:3) and repeat the test. The result must be multiplied by the dilution used.

11. LIMITATIONS OF THE PROCEDURE

A serum sample obtained during the acute phase of infection, when only IgM antibodies are present, may be negative by this procedure.

The Rubella IgM level should be determined using Rubella IgM kit (REF 81031). Alternatively, a second serum sample obtained 8-14 days later, should be tested in parallel to determine an increase in the IgG antibody level.

The test result should be used in conjunction with information available from the evaluation of the case history or other diagnostic procedures.

Although the control serum contained in the kit is calibrated according to the WHO Reference Serum, some discordance may be observed when the serum is tested with different analytical techniques.

12. ANALYTICAL SPECIFICITY

The presence in the test sample of IgG antibodies against viruses such as Cytomegalovirus, Epstein Barr., Herpes Simplex, Measles and Mumps does not influence the result of the test

13. DIAGNOSTIC SENSITIVITY AND SPECIFICITY

242 samples were tested in an experimentation: 25 were negative, 217 positive. The samples were tested with another commercial method. The results are summarized in the following table.

Table 1.

		REFERENCE	
		+	-
DIESSE	+	216	2
	-	1	23

The kit shows a specificity of 92.0% (CI 95%: 75.0–97.8) and a sensitivity of 99.5% (CI 95%: 97.4–99.9).

14. PRECISION

Table 2. Within-in run precision

Sample in triplicate	Mean IU/ml	SD	CV%
Serum 1	5.9	0.21	4
Serum 2	5.3	0.23	4
Serum 3	77.5	3.75	5
Serum 4	30.8	2.75	9
Serum 5	17.3	2.67	15
Serum 6	28.3	2.58	9

Table 3. Between-run precision

	IU/ml						Mean	SD	CV%
	Run 1	Run 2	Run 3	Run 4	Run 5	Run 6			
Serum 1	5.9	6.4	6.2	6.5	5.0	5.0	5.8	0.68	12
Serum 2	5.3	5.0	5.2	5.2	5.0	5.0	5.1	0.13	3
Serum 3	77.5	74.8	77.9	78.2	73.6	71.7	75.6	2.66	4
Serum 4	30.8	30.8	29.4	31.6	26.7	26.1	29.2	2.31	8
Serum 5	17.3	18.2	15.4	14.7	12.1	13.5	15.2	2.29	15
Serum 6	28.3	32.6	31.7	32.2	27.0	30.0	30.3	2.27	7

15. REFERENCES

1. G.B. Wisdom; Enzyme-Immunoassay. Clin. Chem. 22: 1243 (1976)
2. M. Gravell, P. Dorsett et al; Detection of antibody to rubella virus by enzyme-linked immunosorbent assay. J. Inf.Dis. 136 (suppl.) S300 (1977)
3. A. Voller and D. Bidwell; A simple method for detecting antibodies to rubella. Brit. J. Exp. path. 56: 338 (1975)

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